

RHAPSODY RECUMBENT BATHING SYSTEM



ArjoHuntleigh bathing systems ensure an optimum environment to support high quality care and good working conditions throughout the assisted bathing cycle.

Our options enhance the bathing experience with the deep relaxing effects of *Hydromassage*, the sensitive cleaning action of *Hydrosound* and the gentle sensory stimulation of *Sound & Vision*.

Having an effective way to combat

the ever-present risk of cross-infection is vitally important. Our bathing systems offer easy-to-clean surfaces and an effective, well-proven integrated disinfection system.

Hygiene lifters are an integral part of the bathing system, providing safe and comfortable support for the patient during bathing and transport. ArjoHuntleigh offers four hygiene lifters to suit different needs and requirements.

...with people in mind

www.ArjoHuntleigh.com


KEY FEATURES



Sound & Vision™
Colourful underwater lighting and soft music stimulus enhances the power of bathing to calm stressed or anxious patients and residents.



Hyromassage™
A built-in *Hydromassage* unit creates millions of small bubbles and agitation in the bathwater.

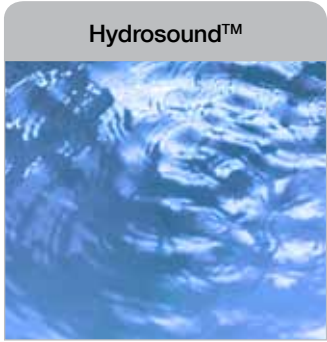

Hydrosound™
Hydrosound offers a gentle touch for cleaning sensitive skin.



P300
The panel that offers the most advanced comfort and safety features:

- Auto-fill+, three preset levels with possible afterfilling
- Three-sensor scalding protection
- Pre-programmed liquid dispensing system
- Digital adjustment of water temperature
- Digital adjustment of *Hydrosound*
- Digital temperature display
- Disinfectant container warning indicates when disinfectant is used up



Product information		
Rhapsody Short	Standard	Long
Total weight, kg (lbs) 390 (860)	470 (1036)	525 (1157)
Water consumption*, liter 185	260	310
Filling & Emptying time at 300 kPa dyn. pressure (Hw 80°C)*, min. 3.00	4.50	5.15
Emptying time, min. 3.00	4.50	5.25
Total length 1970mm (77 ½')	2270mm (89 ½')	2520mm (99 ¼')

Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice.
© and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2012.